

SCHILLER

CARDIOVIT FT-2

SCHILLER's answer to the highest hospital demands

001
Smith John | 21 years 4 months | Male



60



25 mm/s, 10 mm/mV
ECG quality good
An ECG recording can now be taken

Patient

Manual




SCHILLER
The Art of Diagnostics

CARDIOVIT FT-2

Designed for the hospital:

- Advanced cybersecurity features
- Extensive connectivity
- Optional ultrasound spirometry and A4 thermal printer

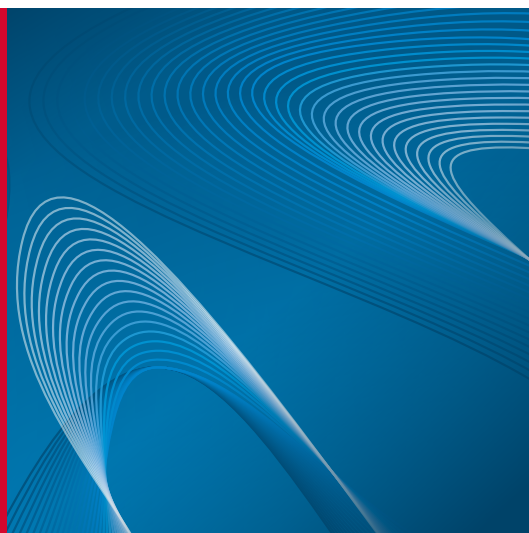

DESIGNED FOR HOSPITAL USE

HIGHLY HYGIENIC: The seamless surface of the large touchscreen is easy to clean. No keyboard is necessary to operate the device. All cables are neatly stowed away in the trolley.

ADVANCED CYBERSECURITY:

- Security-hardened Linux kernel to minimize the risk in case of cyberattacks
- Configurable access control with user name, password and privileges to prevent unauthorized access
- Wi-Fi protocols including certificate-based authentication to ensure secure access to the hospital IT infrastructure

EXTENSIVE CONNECTIVITY:

Bidirectional communication for easy data access via worklist or PDQ, and fast transmission of all reports to EMR/HIS systems.

CUSTOMIZABLE: The CARDIOVIT FT-2 provides an optional thermal printer, allowing to choose the best configuration for paper-based or paperless hospital workflow.

HIGH-END DIAGNOSTIC TOOLS

- Resting ECG with the SCHILLER ECG interpretation algorithm ETM including ETM Sport, automatic interpretation for athletes
- Resting rhythm recording up to 20 minutes
- Spirometry, based on the ultrasound technology of the SpiroScout SP plus (optional)

ADD SPIROMETRY TO YOUR CARDIOVIT FT-2

With the SpiroScout SP plus (optional), the CARDIOVIT FT-2 also becomes a spirometer. Thanks to GANSHORN's ultrasound technology, spirometry recordings are highly accurate and the sensor is maintenance- and calibration-free.

The spirometry option offers FVC, SVC and MVV with a large selection of reference standards, pre and post tests, as well as interpretation of the FVC results. The device provides feedback on the procedure's quality during the test and after each trial. An incentive screen is included for trials with children.



Ultrasound spirometer
SpiroScout SP plus

Large 15" multi-touch display: highly hygienic, easy to clean and easy to use




Anatomical hook up adviser, colour-coded lead quality check and lead reversal detection facilitate ECG acquisition

Spirometry compliant with 2019 ATS/ERS guidelines

Large surface for writing



Optional high resolution A4 thermal printer

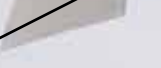

Cables tidily stowed away



Dedicated basket for accessories storage

Easy-glide wheels for smooth manoeuvring





Americas
SCHILLER Americas Inc.
Doral, Florida 33172

North America:
Phone +1 786 845 06 20
Fax +1 786 845 06 02
sales@schilleramericas.com
www.schilleramericas.com

Latin America & Caribbean
Phone +1 305 591 11 21
Fax +1 786 845 06 02
sales@schilleramericas.com
www.schilleramericas.com



France
SCHILLER Médical
F-67160 Wissembourg
Phone +33 3 88 63 36 00
info@schiller.fr
www.schiller.fr



France (distribution France)
SCHILLER France S.A.S.
F-77608 Bussy St Georges
Phone +33 1 64 66 50 00
contact@schillerfrance.fr
www.schiller-france.com



Serbia
SCHILLER d.o.o.
11010 Beograd
Phone +381 11 39 79 508
info@schiller.rs
www.schiller.rs



Slovenia
SCHILLER d.o.o.
2310 Slovenska Bistrica
Phone +386 2 843 00 56
info@schiller.si
www.schiller.si



Asia
SCHILLER Asia-Pacific / Malaysia
52200 Kuala Lumpur, Malaysia
Phone +603 6272 3033
sales@schiller.com.my
www.schiller-asia.com



Germany
SCHILLER Medizintechnik GmbH
D-85622 Feldkirchen b. München
Phone +49 89 62 99 81 0
info@schillermed.de
www.schillermed.de



Spain
SCHILLER ESPAÑA, S.A.
ES-28232 Las Rozas/Madrid
Phone +34 91 713 01 76
schiller@schiller.es
www.schiller.es



Austria
SCHILLER Handelsgesellschaft m.b.H.
A-4040 Linz
Phone +43 732 70 99 0
sales@schiller.at
www.schiller.at



India
SCHILLER Healthcare India Pvt. Ltd.
Mumbai – 400 059, India
Phone +91 22 6152 3333 / 2920 9141
sales@schillerindia.com
www.schillerindia.com



Switzerland
SCHILLER Schweiz AG
CH-8912 Obfelden
Phone +41 44 744 30 00
sales@schiller-schweiz.ch
www.schiller-schweiz.ch



China
Alfred Schiller (Beijing) Medical Technology
Co., Ltd.
100102 Beijing, China
Phone +86 10 84565453
info@schillerchina.com
www.schillermedical.cn



Poland
SCHILLER Poland Sp. z o.o.
PL-02-729 Warszawa
Phone +48 22 843 20 89 / +48 22 647 35 90
schiller@schiller.pl
www.schiller.pl



Turkey
SCHILLER TÜRKIYE
Okmeydani-Sisli – Istanbul
Phone +90 212 210 8681 (pbx)
Fax +90 212 210 8684
info@schiller.com.tr
www.schiller-turkiye.com



Croatia
SCHILLER d.o.o.
10000 Zagreb
Phone +385 1 309 66 59
info@schiller.hr
www.schiller.hr



Russia & C.I.S.
AO SCHILLER.RU
119049 Moscow, Russia
Phone +7 (495) 970 11 33
mail@schiller.ru
www.schiller.ru



United Kingdom
SCHILLER UK Ltd.
Anstruther, Fife KY10 3H
Phone +44 16 190 72 565
sales@schilleruk.com
www.schilleruk.com

SCHILLER reserves the right to make changes in design, specifications and features shown or discontinue the product described at any time without notice or obligation. Images may be representative. Some products may not be regulatory cleared or released for sale in all markets. Contact your SCHILLER representative for the most current information and availability.

